


## Innovative Wound Care Solutions

At Histologics, we are passionate about improving wound care. Our portfolio of wound care devices is designed to aid in treatment of wounds by precision, near-painless debridement, tissue collection for biopsy, and preparation of the wound-bed for advanced treatments. Our devices are FDA registered as brush-curettes for surgical and mechanical debridement and as tissue containers with a snap-off head that can be placed in a vial and shipped for lab analysis of entrapped tissue and fluids. Visit our website for complete details and instructions for use.

### Soft K-Bride®

#### Designed For:

- Wound Bed Preparation
- Small to Large Sized Wound Debridement
- Safety and Precision to Avoid Removing Healed Base
- Tissue Sampling and Biopsy
- Wound Measurement & Sizing

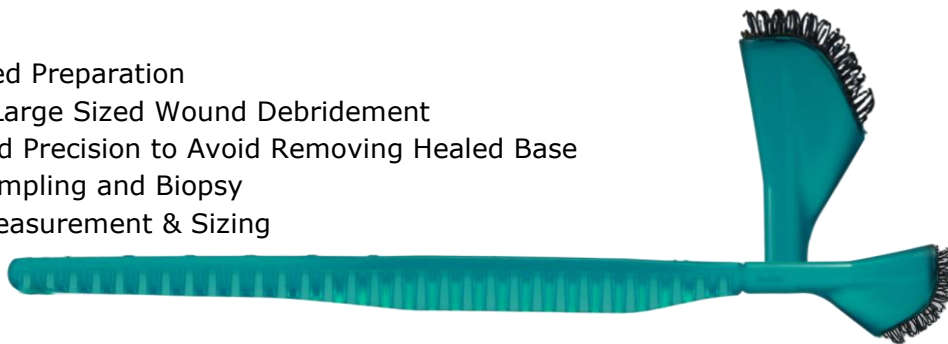

### Soft K-Cot®

#### Designed For:

- Wound Bed Preparation
- Medium to Large Sized Wound Debridement
- Tactile, Safe Wound Debridement



### Soft K-Rette®

#### Designed For:


- Crevice Wound Debridement, Visible Tunnels
- Curettage
- Tissue Sampling and Biopsy




### Soft Biopsy®+D

#### Designed For:

- Wound Tissue Sampling and Biopsy
- Small Wound Debridement
- Small Wound-Bed Preparation

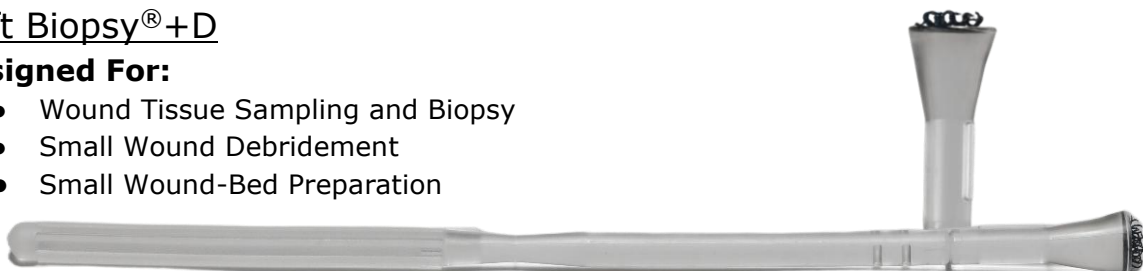

Histologics wound care devices feature our patented Kylon® medical fabric. Kylon® is a nylon “hooked” fabric configured as an array of micro-curettes that can gently cut into soft tissue allowing millimeter-by-millimeter depth precision. This enables care-givers to remove slough and biofilm while stimulating the wound without removing viable tissue. A key benefit is the ability to

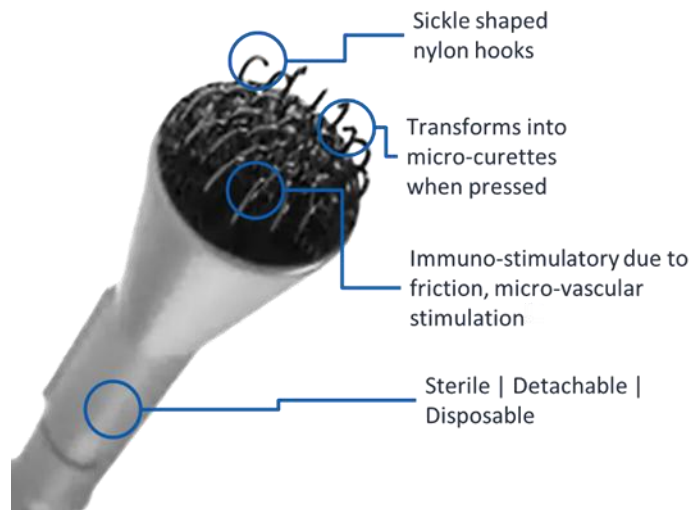

impart agitation to underlying tissue through micro-punctation and frictional heat which is intended to impart blood flow and immune stimulation.

Kylon® devices enable wound-bed preparation and debriding with precise and gentle control avoiding trauma that is common with scalpel or traditional sharp curettes.

## Product Features

- Simple and intuitive design and usage
- A Snap-off head enabling easy transport to the Lab
- 3 year shelf life
- Measurement increments for wound sizing (Soft K-Bride®)
- Sterile disposable plastic, individually packaged
- FDA Compliant and Registered
- Compatible reimbursable procedures

## CPT Codes for Procedure Reimbursement

- CPT 97602 - Non-selective debridement
- CPT 97597/8 - Selective debridement
- CPT 11042/5 - Debridement, subcutaneous tissue
- CPT 11043/6 - Debridement, muscle/fascia
- CPT 15002/3 - Surgical prep for skin graft

## Website



## Order form

